

# CONFERENCE™

A conversation amplifier with two microphones and an inductive loop built into an elegant folder.

## Practical conference folder

The Conference™ is a unique product which enables people with impaired hearing to participate in meetings at work or at home. Two microphones, an amplifier and an inductive loop are built into an elegant folder, which the user places on the table in front of him/her.

## Dual microphones

Conference is equipped with two microphones that shorten the distance to the meeting attendants and transmit the sound to the hearing instrument. If you want even shorter distance to the person speaking, you can use an external microphone.



## Built-in loop

The built-in inductive loop makes it possible for hearing instrument users to use the T mode on the hearing instrument. A neck loop is also supplied, which gives more freedom of movement. The internal loop automatically disconnects when the neck loop is connected.


## Record a meeting


Sometimes, when not wanting to miss out on important information, it is practical to record the conversation from a meeting. The Conference is equipped with a 3.5 mm jack for both recording and playback in mono.

## Runs on batteries

The Conference runs on batteries and has low power consumption. Simply change batteries when they run out.

# TECHNICAL DATA

Frequency range:	200-8000 Hz, -3 dB rel 1 kHz.
Microphones:	2 pc built-in.
Inputs:	External microphone.
Outputs:	Neck loop. Tape recorder (3.5 mm, input: tip, output: ring).
Power supply:	4 pcs LR6/AA 1.5 V penlight.
Battery life:	Approx. 100 hours.
Field strength:	+12 dB rel 100 mA/m at 0.7 mV (70 dB SPL) from the microphones at a distance of 45 cm from the loop.
Distortion:	<3%.
Gain:	0.7 mV from the microphones produces 250 mV out at the tape recording. The tape contact can also be used for playback.
Current consumption:	25-200 mA.
Weight:	Approx. 1.3 kg.
	Any issues relating to the EU Medical Device Directive 93/42/EEC should be directed to GN ReSound A/S.
Ordering number:	415 04 04

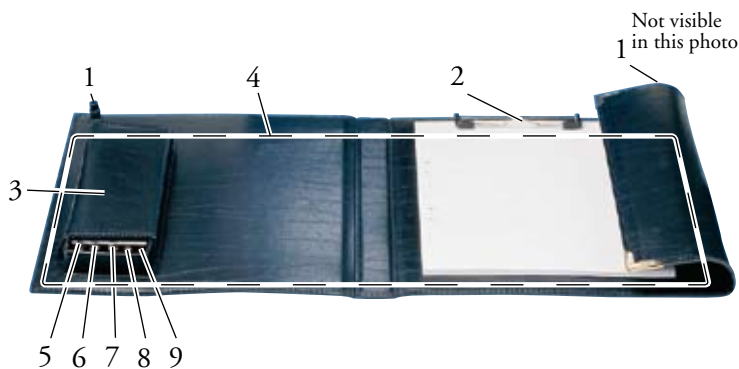

1. Microphones
2. Paper clip
3. Battery and neck loop compartment
4. Built-in inductive loop
5. Tape recorder input/output
6. Microphone input
7. Neck loop output
8. ON/OFF switch
9. ON/OFF LED indicator

#### CONTACT:

**GN ReSound AB ALD Division**  
Svarvarg. 9, 2 tr  
SE-112 49 Stockholm  
Sweden

Tel. +46 8 441 16 70  
Fax +46 8 441 16 80  
info.ald@gnresound.se  
www.gnresound-ald.com

#### MANUFACTURER

*in the sense of EU Medical  
Device Directive 93/42/EEC:*  
**GN ReSound A/S**  
Lautrupbjerg 7  
DK-2750 Ballerup  
Denmark

*Service and product-related  
questions: see "contact".*