

t-hook

A portable handsfree for hearing instrument users

Handsfree for hearing instrument users

T-hook is a light and comfortable handsfree for hearing instrument users who are using mobile (cellular) phones or cordless telephones (i.e. DECT).

T-hook can be used with both in-the-ear (ITE) and behind-the-ear (BTE) hearing instruments. The design makes T-hook easy to place behind the ear and it is comfortable to wear. The sound is transferred to the hearing aid using the T-switch mode.

With T-hook you experience excellent sound quality. T-hook is easy to use and has a microphone incorporated in the handsfree cable.




Fits most mobile and cordless phones

T-hook has a 2.5 mm plug and can, with adapters, be used with almost all mobile and cordless telephones. Adapters are available from GN ReSound.


Some mobile phones and cordless phones are equipped with a 2.5 mm connector, eliminating the need for an adapter.




HOW TO USE THE T-HOOK

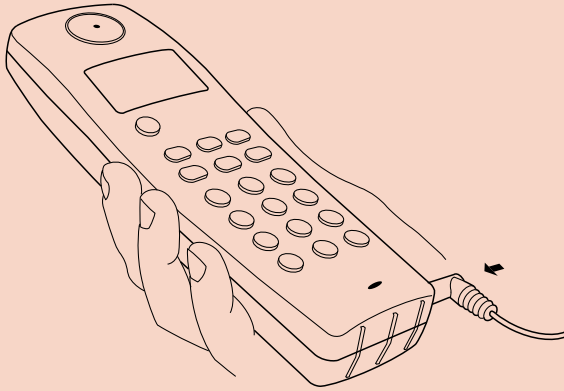

Direct connection to the telephone.

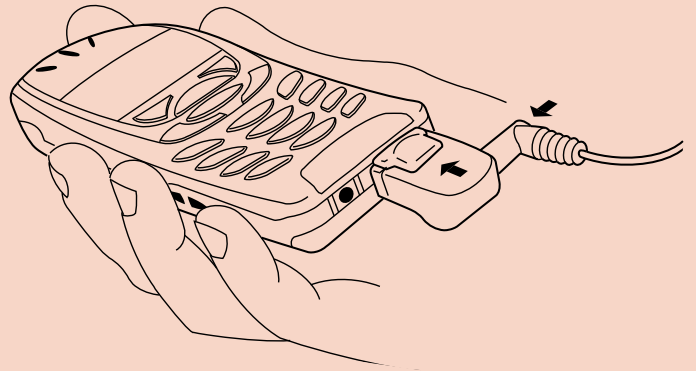

Connection by means of an adapter.

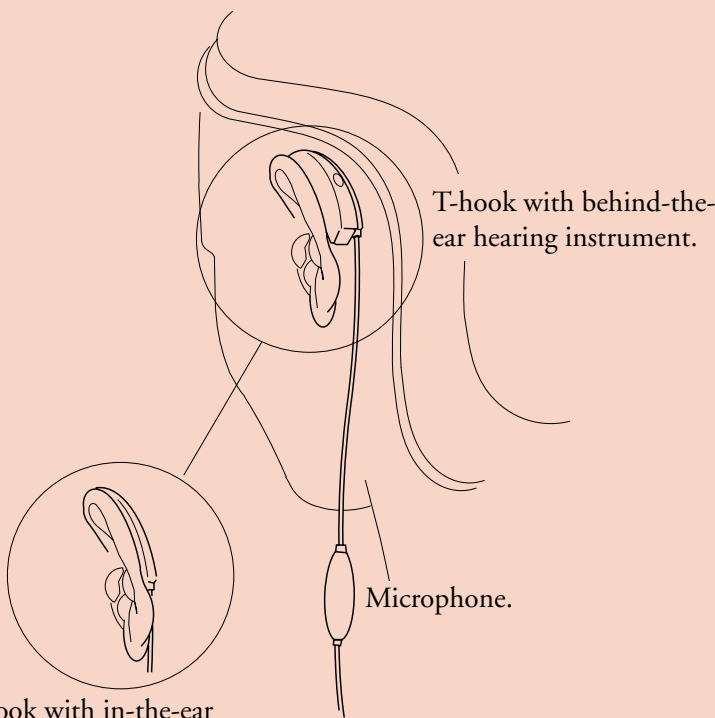

T-hook with behind-the-ear hearing instrument.

T-hook with in-the-ear hearing instrument.

Technical data

Impedance:	177 Ohm.
Resistance:	175 Ohm.
Inductance:	4.4 mH.
Cut-off frequency, -3 dB:	6 kHz.
Sensitivity:	0.6 V in gives a field strength of 100 mA/m at a distance of 5 cm.

CE Any issues relating to the EU Medical Device Directive 93/42/EEC should be directed to GN ReSound A/S.

Ordering number "T-hook Universal 2.5 mm": 200 21 37



T-hook has received the award "Excellent Swedish Design".

CONTACT:

GN Resound AB ALD Division
 Svarvarg. 9, 2 tr
 SE-112 49 Stockholm
 Sweden
 Tel.: +46 8 441 16 70
 Fax: +46 8 441 16 80
 info.ald@gnresound.se
 www.gnresound-ald.com

MANUFACTURER

in the sense of EU Medical Device Directive 93/42/EEC:
GN Resound A/S
 Lautrupbjerg 7
 DK-2750 Ballerup
 Denmark
 Service and product related questions: contact GN ReSound ALD Division.